

# Podoblock™

Innovative vet & solutions  
through experience

Catalogue

2018





10:30  
Call Joachim,  
feedback on  
MagRayCloxx





# Podoblock B.V.

## We are Podoblock!

...and pleased to offer you a new look into our comprehensive product line. All our products have been developed in order to help you make quick and safe, high-quality X-rays. Every piece of equipment in this catalogue has been designed, built and tested by us, based on experience, earned during 25 years of working professionally with horses. All made in the Netherlands.

They are made of the highest quality materials, built to last for decades. Podoblock supplies accessories for every X-ray system on the market, both for standard and special tailor-made configurations.



We supply our products through a network of distributors or, if needed, directly to you. Through solid connections with our distributors and vets around the world we built a network for coming up with the best ideas and the latest problems to solve.

Being the no-nonsense company that we always have been, we continue to take pride in knowing our customers personally and use short communication lines.

'Impossible' is an adjective that we abhor and we will release all our intellectual and critical faculties in order to solve any problem you run into. When you contact us you will talk to one of us - and no one else - and you will deem us honest, trustworthy and reliable.

Feel free to roam around and don't hesitate to contact us if you have any questions.

Any!

*Alex Berends*

### Company address:

Verlengde Nijverheidsweg 10  
9482 WT Tynaarlo  
The Netherlands

☎ Alex  
Office  
Fax

0031-651 828 117  
0031-508 200 257  
0031-508 200 212

Web: [www.podoblock.com](http://www.podoblock.com)  
@: [info@podoblock.com](mailto:info@podoblock.com)

Chamber of Commerce: 58189998  
VAT nr: NL852918227B01



# Contents

		Page
H1	 Podoblock	6
H2	 CR Cassette Solutions	10
H3	 Dental X-ray Solutions	15
H4	 DR Flatpanel Solutions	19
H5	 Flatpanel Manufacturers	22
H6	 Ultrasound Solutions	29
H7	 Generator Solutions	34
H8	 Protection & Safety	40
H9	 X-ray Markers	46
H10	 Special X-ray Solutions	48
H11	 Miscellaneous	50




At Podoblock we recognize the challenges that equine vets meet on a day-to-day basis. From the onset of our business in 2008 we have been aware of the hurdles that come with equine X-raying.

Our very first contraption, the Podoblock, has been a whopping success ever since. It always guarantees the same angle in every position and it minimizes radiation in a maximal fashion, because of the hands-free system.

The Podoblock is robust, solid and almost indestructible. It comes in 3 different heights and has a number of optional extras and add-ons, laid out for you on the following pages.

# Podoblock

- Made of aluminium
  - Weight Podoblock: 6,5 kg
  - 4 rubber feet
  - includes 1 Topblock of your choice
  - LM / AP / 45° Handsfree system
- 
- To be used in combination with our Podoholder, Podo-Protectionbox or Cowboybox, for optimal radiation distance
  - Avoids fingers on your X-ray
  - Height from ground to tip-off hoof in Podoblock is either 90, 120 or 170 mm

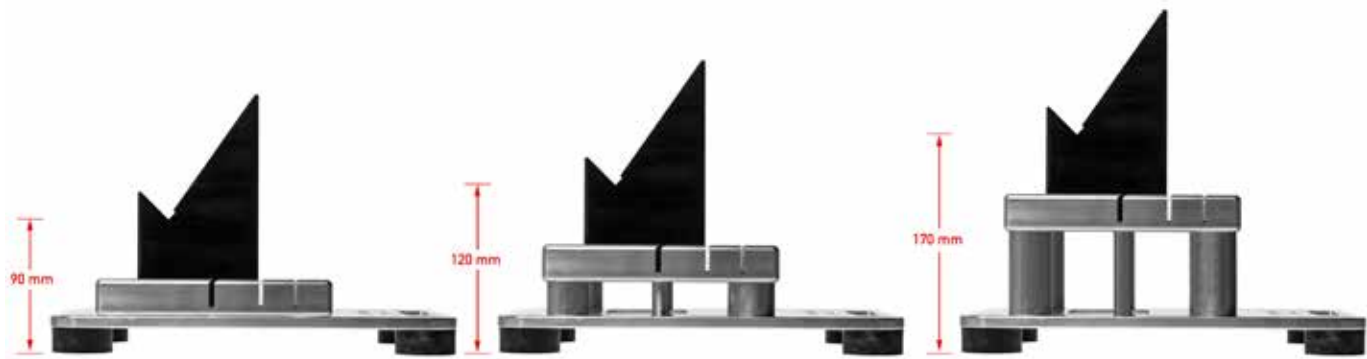

Product	Ordernumber	Price
Podoblock 90 mm	PB 01AL090	€ 475,-
Podoblock 120 mm	PB 01AL120	€ 475,-
Podoblock 170 mm	PB 01AL170	€ 475,-

## Topblock + optional Wedge

Easily interchangeable blocks for our Podoblock  
Angle diversity 45° to 60° with 5° increments

- Made of synthetic PE500
- Available in black only

Product	Ordernumber	Price
Topblock 45°	TB 0250045	€ 150,-
Topblock 50°	TB 0250050	€ 150,-
Topblock 55°	TB 0250055	€ 150,-
Topblock 60°	TB 0250060	€ 150,-
Topblock 60° Upright Pedal	TBD 0250060	€ 150,-
Topblock Wedge 10°	TBW 0250010	€ 100,-



## Topblock Flatplate

For LM / AP / 45° takes with the foot standing flat or upright pedal

- Made of synthetic PE500
- Available in black only



Product	Ordernumber	Price
Topblock Flatplate AP	TBF 02500AP	€ 150,-
Topblock Flatplate LM / AP	TBF 02500LMAP	€ 150,-
Topblock Flatplate LM / AP 45°	TBF 0250045	€ 150,-

## Optional Bottomplate

Replaces the rubber feet of the regular Podoblock in order to create sliding effect

- Made of HMPE500
- Weight: 500 gr
- Thickness of the plate is 10 mm



Product	Ordernumber	Price
Bottomplate	PV 01PE500	€ 75,-

## MagRayCloxx

Magnetic X-ray block that greatly enhances both your and the horse's mobility and working speed. The blocks stay on the horse while making oblique X-rays.

- Made of synthetic PE500
- Clean break-over on front and side
- Magnets to fix to the horseshoe
- Comes with radiolucent foil for extra fixation



Design:  
Dr. Joachim Onderka



Product	Ordernumber	Price
2 Blocks / 1 Hand dispenser / 5 Foil rolls	MRC 01SET2-5	€ 395,-
4 Blocks / 1 Hand dispenser / 12 Foil rolls	MRC 01SET4-12	€ 725,-
1 Hand dispenser / 12 foil rolls	MRC 01FOIL1-12	€ 65,-



## Angleator

A nifty and extremely handy device for X-rays of the equine foot

With a stable base in combination with a rotating disk, this contraption will facilitate all your LM / AP hoof exams. There is no more need to fidget and fumble with the horse's foot, foot-block or X-ray source

The Set includes a dummy for the other foot. Base and dummy can also be ordered separately

The Angleator comes in a handy cotton bag for transportation

The plateau has predefined angle-markings for your convenience and has pre-incised slots at 0°, 45° and 90°. The red and green crosslines are accommodated with lead markers

- Made of HMPE500 and ABS
- Height: 82 mm
- Weight of set: 5,0 kg
- Diameter plateau: 20 cm
- Total Width: 27,5 cm
- Total Length: 34 cm
- Anti-slip surface



Product	Ordernumber	Price
Angleator Set	AK 10ADSET	€ 800,-
Base	AK 10ADBASE	€ 620,-
Dummy	AK 10ADDUM	€ 180,-



The Angleator is best used in combination with our Podoholder, Podo-Protectionbox and Cowboybox. These items are to be ordered separately.



# 2

## CR Cassette Solutions



Computed radiography (CR) systems are still much in use in the veterinary world. At Podoblock we offer a complete product line in order to enable your CR system for your equine tasks, outlined for you on the following pages.

## Podoholder

To be used in combination with our Podoblock and Combibox  
Protects your cassette and always secures the right positions  
for standard hoof exams

Optional to the Podoholder are the Podoholder Lead and  
Podoholder Grid. If you want to use a Grid or Backscatter lead,  
please contact us for optional possibilities



Size	Ordernumber	Price
18 x 24	PH 03182418	€ 110,-
18 x 24 Horizontal	PHB 03182418	€ 110,-
24 x 30	PH 03243018	€ 110,-
24 x 30 Horizontal	PHB 03243018	€ 110,-
35 x 43 Horizontal	PHB 03354318	€ 110,-
8" x 10"	PH 03232818	€ 110,-
8" x 10" Horizontal	PHB 03232818	€ 110,-
10" x 12"	PH 03283318	€ 110,-
10" x 12" Horizontal	PHB 03283318	€ 110,-

For cassette measurements please check table on p 14

## Knee Holder

Easy to use distance holder - especially for knee exams -  
that creates maximum distance and safety

Standard supplied with a telescopic stick that can be  
adjusted to the left or the right in a split second

- Made of 3 mm polycarbonate.
- For use with all Analogue and CR cassettes



Size	Ordernumber	Price
18 x 24	KH 04182418	€ 250,-
24 x 30	KH 04243018	€ 250,-
8" x 10"	KH 04232818	€ 250,-
10" x 12"	KH 04283318	€ 250,-

For cassette measurements please check table on p 14





## Distance Cassette Holder

Crafted to create maximum distance from the radiation as well as from the horse

Comes with a 3D-adjustable telescopic stick with only one knob for securing



Telescopic stick

- Weight: 1,5 kg
- Made of super strong durable lightweight polycarbonate

Size	Ordernumber	Price
DCH 18 x 24	AC 04182418	€ 300,-
DCH 18 x 24 Fetlock	ACK 04182418	€ 300,-
DCH 24 x 30	AC 04243018	€ 300,-
DCH 35 x 43 Horizontal	ACB 04354318	€ 300,-
DCH 8" x 10"	AC 04232818	€ 300,-
DCH 8" x 10" Fetlock	ACK 04232818	€ 300,-
DCH 10" x 12"	AC 04283318	€ 300,-

For cassette measurements please check table on p 14

## Distance Cassette Holder with foot option

Same as above, but with an easy-click-on and in height adjustable telescopic foot

Telescopic foot available in 2 lengths:

- Short: 28 - 49 cm for hock
- Long: 100 - 190 cm for back/neck
- Weight: 1,9 kg
- Made of super strong, durable, lightweight polycarbonate



Telescopic foot

Size	Ordernumber	Price
DCH 18 x 24 + Short Foot	ACS 04182418	€ 375,-
DCH 24 x 30 + Short Foot	ACS 04243018	€ 375,-
DCH 35 x 43 + Long Foot	ACBS 04354318	€ 375,-
DCH 8" x 10" + Short Foot	ACS 04232818	€ 375,-
DCH 10" x 12" + Short Foot	ACS 04283318	€ 375,-
Optional Support Foot (S)	LS 10300MM	€ 75,-
Optional Support Foot (L)	LS 101200MM	€ 75,-

For cassette measurements please check table on p 14



## Skylinebox

This is the most effective product on the market for high-quality skyline X-rays and makes taking skyline exams simple because of the special built-in angle

After exposure it is extremely easy to open and to reload the tunnel without replacing the horse



The Skyline tunnel is available for all sizes of analogue and CR cassettes

- Made of polycarbonate 10 mm
- 4 rubber feet and an anti-slip surface
- Skyline angle 5°

Size	Ordernumber	Price
18 x 24	SK 05182420	€ 275,-
24 x 30	SK 05243020	€ 275,-
8" x 10"	SK 05232820	€ 275,-
10" x 12"	SK 05283320	€ 275,-

For cassette measurements please check table on p 14

## Combibox

The Combibox has all the features of the Skyline tunnel above and can additionally be used for LM / AP flat and upright pedal takes, when turned over



With the optional Podoholder the Combibox creates a handsfree system that allows for making Skyline, LM and AP X-rays of the foot.

For optimal results use it in combination with the Dummy to have both feet at the same height

The Combibox is available for all sizes cassettes, both analogue and CR

- made of polycarbonate 10 and 18 mm
- 70 mm in height
- rubber feet and an anti-slip surface
- Skyline angle 5°




Product	Ordernumber	Price
Combibox 18 x 24 LM	CO 061824081	€ 400,-
Combibox 18 x 24 LM / AP	CO 061824082	€ 450,-
Combibox 18 x 24 LM / AP & Dummy	CO 061824083	€ 625,-
Combibox 24 x 30 LM	CO 062430081	€ 400,-
Combibox 24 x 30 LM / AP	CO 062430082	€ 450,-
Combibox 24 x 30 LM / AP & Dummy	CO 062430083	€ 625,-
Combibox 8" x 10" LM	CO 062328081	€ 400,-
Combibox 8" x 10" LM / AP	CO 062328082	€ 450,-
Combibox 8" x 10" LM / AP & Dummy	CO 062328083	€ 625,-
Combibox 10" x 12" LM	CO 062833081	€ 400,-
Combibox 10" x 12" LM / AP	CO 062833082	€ 450,-
Combibox 10" x 12" LM / AP & Dummy	CO 062833083	€ 625,-
Dummy LM	CO 06DUM081	€ 200,-

For cassette measurements please check table on p 14

## Measurements X-ray cassettes - Analogue & Digital CR

Please check the cassette measurements in case you're not sure what type of cassette you use. We supply cassette holders for every X-ray system. If you have different measurements cassettes, please contact us

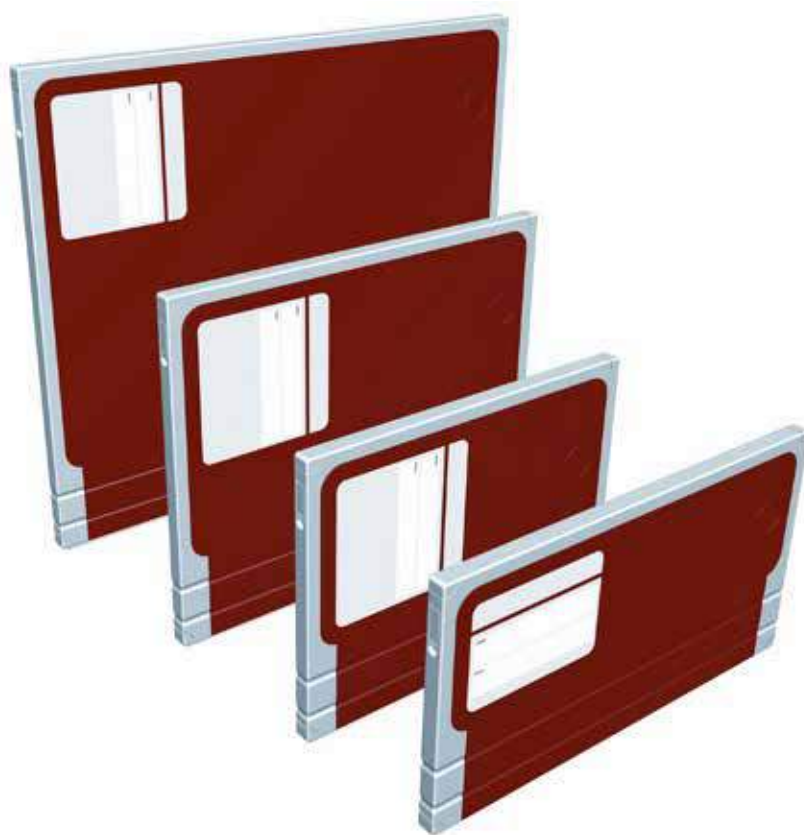

Cassette size (mm)	Film size (inch)	Film size (cm)	Ordernumber
208 x 267		18 x 24	... ..1824...
267 x 327		24 x 30	... ..2430...
385 x 385	14" x 14"	35 x 35	... ..3535...
385 x 460	14" x 17"	35 x 43	... ..3543...
460 x 460	17" x 17"	43 x 43	... ..4343...
230 x 280	8" x 10"	20 x 25	... ..2328...
280 x 330	10" x 12"	25 x 30	... ..2833...



# 3

## Dental X-ray Solutions



Five years ago we could not have envisioned the demand for dental X-ray solutions we get these days. From the onset it soared and we tried to lead the way from then on with practical solutions

At this junction we have a number of standard products to improve and yet simplify these usually difficult tasks. In case you have specific wishes with regard to new products, we are happy to collaborate with you on possible solutions to your needs

## Denfinity

This new light-weight portable and collapsible stainless steel stand offers you the freedom to make high-quality Skull X-rays in hands-free fashion

The Denfinity stand is a fully stainless steel base plate and vertical telescopic pole that guarantees a secure positioning and right height adjustment

Starting with the base element of the Denfinity, the base plate and pole, you can configure with a number of options that can easily and directly be placed on the Denfinity base

The aluminium 3D-joint is used in combination with our Protectionbox

The LM-Head support is best used in combination with our Protectionbox, Podo-Protectionbox and Cowboybox and also fixates your CR cassette or DR flatpanel

The Protectionbox secures your Cassette or Flatpanel when use in combination with the 3D-joint and LM-Head support

When ordering the optional Protectionbox or Podoholder, be sure to state the required cassette measurements you use



## Denfinity options



### Base Element

Base plate: 350 x 405 x 90 mm  
Vertical pole: 960 mm  
Vertical travel: 600 mm  
Material: Stainless steel



### 3D-joint

Aluminium 3D-swivel joint used in combination with Dental- or Protectionbox



### LM-Head support

Polycarbonate 10 mm  
Swivel range: 180°  
Universal adapter for fixation Podoholder or Protectionbox



### Headrest

Polycarbonate 10 mm  
Full radiolucent, no metal parts in detection area



### Protectionbox

Material: Polycarbonate 3 mm  
Fully closed cassette or DR-panel Protectionbox  
*Requires 3D-joint to function*

Product	Ordernumber	Price
Base Element	DSJ 1050401	€ 550,-
3D-Joint	DSJ 10DNL3D	€ 425,-
LM Head-support	DSJ 10LMHEAD	€ 350,-
Headrest	DSJ 10HREST	€ 200,-
Protectionbox	DSJ 10PROT	€ 600,-

## XtandR

A revolutionary set of interchangeable biteplates for equine head radiographs

The XtandR is only intended for taking X-rays of the incisors

Is compatible with all available speculi on the market

- Material: HMPE1000
- Can withstand a pressure of 5,13 KN / 513 kg



Radiolucent bite plates sample Xray



### Note:

The speculum shown is **not** provided with the XtandR

### Attention:

In case you have a CAPPS speculum, be aware there is a separate ordernumber for your XtandR

Product	Ordernumber	Price
XtandR	RBP 10BPPZ	€ 650,-
XtandR (CAPPS model)	RBP 10CAPPS	€ 650,-



## The Juliuster

Designed, tested and approved in collaboration with the European Equine Dental Consultancy (EEDC )

Easy-to-use device to offset the lower jaw to the left or right for open-mouth views with offset-technique

Great to investigate the roots and the occlusion surface of the cheek teeth. Dorso-ventral views can be made with 5 cm offset to the left and right

- Material: HMPE500
- Offset: 5 cm left/right
- Adjusted in 12 mm increments
- Strap: 148 cm Biothane



Product	Ordernumber	Price
Juliuster	DRS 1011560	€ 375,-

## Dental X-Block

Can be used directly between the incisors to quickly create 60, 80 or 100 mm space

- Material: Rubber 95° shore
- Weight: 650 gr
- Size: 60 x 80 x 100 mm



Product	Ordernumber	Price
Dental X-Block	DXB 106810	€ 70,-

# 4

## DR Flatpanel Solutions



Digital Radiography (DR) techniques have evolved rapidly in a relatively short period of time, not just concerning software, but also hardware, maybe even more so. The first DR flatpanels were a heavy, wired and inconvenient version of what they are today, yet suitable for most equine vets, thanks to our products.

The evolution to thin, light-weight, wireless and easily manageable flatpanels was a rapid one and has been a constant challenge for us to keep up with the changes and developments on the market and provide adequate solutions.

With our DR-product package we believe we have the most appropriate solution for today's equine vet. At least we can't envision a practice without it!

## DR-Protectionbox

Light-weight protectionbox - Handheld

Depending on your DR system, the box comes with a mounted grip handle or with our Greppstr 2.1

Optional: - Grid Compartment (tailor-made)



## DR-Podo-Protectionbox

Protectionbox with 5 mm adaptors to place in our Podoblock and Combibox, that guarantees the same direction and angle with hoof examinations at all times

Depending on your DR system, the box comes with a mounted grip handle or with our Greppstr 2.1

Optional: - Grid Compartment (tailor-made)

## DR-Cowboybox

Protectionbox, Podo-protectionbox and Remote cassette holder all-in-one

Standard with mountable telescopic grip, which can be turned in 3D positions to get extra distance and a mountable and height-adjustable 300 mm support foot

Depending on your DR system, the box comes with a mounted grip handle or with our Greppstr 2.1

Optional: - Grid Compartment (tailor-made)



## DR-Skylinebox

Because of the special angle the horse is standing comfortably and there is more space for the generator

Top surface is supplied with radiolucent anti-slip surface and the door of the tunnel is secured by a lock

● Material: 20 mm polycarbonate

## DR-Sleeve

Soft cover Sleeve made of PVC with Velcro closing

Front side is fully radiolucent with marked detection area

Inside foam at the edges and corners for extra protection

Depending on your DR system, the Sleeve comes with a mounted grip handle or with our Greppstr 2.1



Product	Ordernumber	Price
DR-Protectionbox	See Pages 22-28	€ 450,-
DR-Podo-Prot. box	See Pages 22-28	€ 450,-
DR-Cowboybox	See Pages 22-28	€ 750,-
DR-Skylinebox	See Pages 22-28	€ 525,-
DR-Sleeve	See Pages 22-28	€ 400,-





## DR-Combibox LM

The Skyline side of the Combibox has the same features as the original DR-Skylinebox. Because of the special angle the horse is standing comfortably and there is more space for the X-ray generator

When turned over the tunnel has extra features to facilitate LM and AP flat views of the foot

The special adapters on the DR-Combibox allow you to take "hands-free" X-rays in combination with our DR Podo-Protectionbox or DR-Cowboybox

This way the DR flat panel is protected and always in the right position

## DR-Combibox LM / AP

The possibility to choose between only the LM adapter or, your best option of course, both LM and AP adapters

This way it is not necessary to replace the foot for both exams



## DR-Combibox SET

Both standing feet on the same height is advised with LM and AP examinations. That is why our Combibox Dummy can also be supplied with the DR-Combibox

The Combibox Set is included:

- DR-Combibox LM / AP
- Combibox Dummy LM

## DR-Cowboy Package

The optimal combined package of DR Products to make all X-rays needed, both safely and in a standardized fashion

This is what you need, for either mobile or stationary use in your clinic

The package deal includes:

- Podoblock
- DR-Cowboybox
- DR-Combibox LM / AP
- Combibox Dummy LM



Product	Ordernumber	Price
DR-Combibox LM	See Pages 22-28	€ 800,-
DR-Combibox LM/AP	See Pages 22-28	€ 825,-
DR-Combibox SET	See Pages 22-28	€ 990,-
DR-Cowboy Package	See Pages 22-28	€ 2.150,-

# 5

## Flatpanel Manufacturers

### Agfa



**1**

DR10s Wireless



**2**

DR14s Wireless



**3**

DX-D 40C Wireless  
DX-D 40G Wireless



**4**

DX-D 45C Wireless  
DX-D 45G Wireless

Producttype	1	2	3	4
DR-Protectionbox	DR 03ISO2430	DR 03ISO1417	DR 03ISO1417	DRVI 031012N
DR-Skylinebox	DR 05ISO2430	DR 05ISO1417	DR 05ISO1417	DRVI 051012N
DR-Combibox LM	DR 06ISO2430 1	DR 06ISO1417 1	DR 06ISO1417 1	DRVI 061012N 1
DR-Combibox LM / AP	DR 06ISO2430 2	DR 06ISO1417 2	DR 06ISO1417 2	DRVI 061012N 2
DR-Combibox SET	DR 06ISO2430 S	DR 06ISO1417 S	DR 06ISO1417 S	DRVI 061012N S
DR-Podo-Protectionbox	DR 07ISO2430	DR 07ISO1417	DR 07ISO1417	DRVI 071012N
DR-Cowboybox	DR 08ISO2430	DR 08ISO1417	DR 08ISO1417	DRVI 081012N
DR-Cowboy Package	DR 08ISO2430 P	DR 08ISO1417 P	DR 08ISO1417 P	DRVI 081012N P
DR-Sleeve	DR 10ISO2430	DR 10ISO1417	DR 10ISO1417	DRVI 101012N



**1**  
Canon CXDI-31



**2**  
Canon CXDI-60C  
Canon CXDI-60G  
Eklin Mark IIG  
Eklin Mark IIIG  
Scope-XS  
XDR2



**3**  
Canon CXDI-70C  
Canon CXDI-701C  
Canon CXDI-701G  
Agfa DX-D 30C  
Wireless



**4**  
Canon CXDI-80C  
Canon CXDI-801C  
Canon CXDI-801G  
SoundEklin **Sprint** Air  
Agfa DX-D 35C  
Wireless



**5**  
Canon CXDI-710C  
Wireless



**6**  
Canon CXDI-810C  
Wireless

Producttype	1	2	3	4	5	6
DR-Protectionbox	DR31 03323330	DR60 03343824	DR 03ISO1417	DR80 03383015	DR 03ISO1417	DR810 03383015
DR-Skylinebox	DR31 05323330	DR60 05343824	DR 05ISO1417	DR80 05383015	DR 05ISO1417	DR810 05383015
DR-Combibox LM	DR31 06323330 1	DR60 06343824 1	DR 06ISO1417 1	DR80 06383015 1	DR 06ISO1417 1	DR810 06383015 1
DR-Combibox LM / AP	DR31 06323330 2	DR60 06343824 2	DR 06ISO1417 2	DR80 06383015 2	DR 06ISO1417 2	DR810 06383015 2
DR-Combibox SET	DR31 06323330 S	DR60 06343824 S	DR 06ISO1417 S	DR80 06383015 S	DR 06ISO1417 S	DR810 06383015 S
DR-Podo-Protectionbox	DR31 07323330	DR60 07343824	DR 07ISO1417	DR80 07383015	DR 07ISO1417	DR810 07383015
DR-Cowboybox	DR31 08323330	DR60 08343824	DR 08ISO1417	DR80 08383015	DR 08ISO1417	DR810 08383015
DR-Cowboy Package	DR31 08323330 P	DR60 08343824 P	DR 08ISO1417 P	DR80 08383015 P	DR 08ISO1417 P	DR810 08383015 P
DR-Sleeve	DR31 10323330	DR60 10343824	DR 10ISO1417	DR80 10383015	DR 10ISO1417	DR810 10383015

## Fuji FDR



**1**

D-EVO C24i Wireless



**2**

D-EVO II C24 Wireless



**3**

D-EVO G35 Wireless



**4**

D-EVO C35i Wireless



**5**

D-EVO II C35 Wireless



**6**

D-EVO II C43 Wireless

Producttype	1	2	3	4	5	6
DR-Protectionbox	DRFU 03C24i	DRFU 03IIC24	DRFU 03DEVO	DR 03ISO1417	DRFU 03IIC35	DR 03ISO1717
DR-Skylinebox	DR 05ISO2430	DR 05ISO2430	DRFU 05DEVO	DR 05ISO1417	DR 05ISO1417	DR 05ISO1717
DR-Combibox LM	DR 06ISO2430 1	DR 06ISO2430 1	DRFU 06DEVO 1	DR 06ISO1417 1	DR 06ISO1417 1	DR 06ISO1717 1
DR-Combibox LM / AP	DR 06ISO2430 2	DR 06ISO2430 2	DRFU 06DEVO 2	DR 06ISO1417 2	DR 06ISO1417 2	DR 06ISO1717 2
DR-Combibox SET	DR 06ISO2430 S	DR 06ISO2430 S	DRFU 06DEVO S	DR 06ISO1417 S	DR 06ISO1417 S	DR 06ISO1717 S
DR-Podo-Protectionbox	DRFU 07C24i	DRFU 07IIC24	DRFU 07DEVO	DR 07ISO1417	DRFU 07IIC35	DR 07ISO1717
DR-Cowboybox	DRFU 08C24i	DRFU 08IIC24	DRFU 08DEVO	DR 08ISO1417	DRFU 08IIC35	DR 08ISO1717
DR-Cowboy Package	DRFU 08C24i P	DRFU 08IIC24 P	DRFU 08DEVO P	DR 08ISO1417 P	DRFU 08IIC35 P	DR 08ISO1717 P
DR-Sleeve	DR 10ISO2430	DR 10ISO2430	DRFU 10DEVO	DR 10ISO1417	DR 10ISO1417	DR 10ISO1717



## Rayence / Vatech



**1**

Vatech 1210P  
Samsung 1210P  
Agfa DX-D 25C  
XDR2-K, XDR2-M, XDR2-S  
meX+ 1210P



**2**

1012WCA Wireless  
SoundEklon Sprint Bolt  
Sprint™ Ultralight DR  
XDR3  
meX+ 1012WCA



**3**

Vatech 1012WCC / WGC  
Wireless



**4**

1417PCA / PGA Wireless  
1417WCA / WGA Wireless



**5**

1417WGC / WCC Wireless

## O&R XenOR



**6**

XenOR 24CW Wireless

Producttype	1	2	3	4	5	6
DR-Protectionbox	DRSX 03514528	DRSX 031012WCA	DRSX 031012WCC	DRSX 031417PGA	DRSX 031417WGC	DR 03ISO1012
DR-Skylinebox	DRSX 05514528	DRSX 051012WCA	DR 05ISO2430	DRSX 051417PGA	DR 05ISO1417	DR 05ISO1012
DR-Combibox LM	DRSX 06514528 1	DRSX 061012WCA 1	DR 06ISO2430 1	DRSX 061417PGA 1	DR 06ISO1417 1	DR 06ISO1012 1
DR-Combibox LM / AP	DRSX 06514528 2	DRSX 061012WCA 2	DR 06ISO2430 2	DRSX 061417PGA 2	DR 06ISO1417 2	DR 06ISO1012 2
DR-Combibox SET	DRSX 06514528 S	DRSX 061012WCA S	DR 06ISO2430 S	DRSX 061417PGA S	DR 06ISO1417 S	DR 06ISO1012 S
DR-Podo-Protectionbox	DRSX 07514528	DRSX 071012WCA	DRSX 071012WCC	DRSX 071417PGA	DRSX 071417WGC	DR 07ISO1012
DR-Cowboybox	DRSX 08514528	DRSX 081012WCA	DRSX 081012WCC	DRSX 081417PGA	DRSX 081417WGC	DR 08ISO1012
DR-Cowboy Package	DRSX 08514528 P	DRSX 081012WCA P	DRSX 081012WCC P	DRSX 081417PGA P	DRSX 081417WGC P	DR 08ISO1012 P
DR-Sleeve	DRSX 10514528	DRSX 101012WCA	DR 10ISO2430	DRSX 101417PGA	DR 10ISO1417	DR 10ISO1012

## Konica Minolta



**1**

AeroDR 1012HQ  
Wireless



**2**

AeroDR 1417S  
AeroDR 1417HQ  
Wireless

## PerkinElmer



**3**

XRPad 3025™ Wireless



**4**

XRPad 4336™ Wireless

## Vieworks

**5**

VIVIX-1012N Wireless



**6**

VIVIX-1417W Wireless  
VIVIX-1417N Wireless



Producttype	1	2	3	4	5	6
DR-Protectionbox	DR 03ISO1012	DR 03ISO1417	DR 03ISO1012	DRPE 03XR4336	DRVI 031012N	DR 03ISO1417
DR-Skylinebox	DR 05ISO1012	DR 05ISO1417	DR 05ISO1012	DR 05ISO1417	DRVI 051012N	DR 05ISO1417
DR-Combibox LM	DR 06ISO1012 1	DR 06ISO1417 1	DR 06ISO1012 1	DR 06ISO1417 1	DRVI 061012N 1	DR 06ISO1417 1
DR-Combibox LM / AP	DR 06ISO1012 2	DR 06ISO1417 2	DR 06ISO1012 2	DR 06ISO1417 2	DRVI 061012N 2	DR 06ISO1417 2
DR-Combibox SET	DR 06ISO1012 S	DR 06ISO1417 S	DR 06ISO1012 S	DR 06ISO1417 S	DRVI 061012N S	DR 06ISO1417 S
DR-Podo-Protectionbox	DR 07ISO1012	DR 07ISO1417	DR 07ISO1012	DRPE 07XR4336	DRVI 071012N	DR 07ISO1417
DR-Cowboybox	DR 08ISO1012	DR 08ISO1417	DR 08ISO1012	DRPE 08XR4336	DRVI 081012N	DR 08ISO1417
DR-Cowboy Package	DR 08ISO1012 P	DR 08ISO1417 P	DR 08ISO1012 P	DRPE 08XR4336 P	DRVI 081012N P	DR 08ISO1417 P
DR-Sleeve	DR 10ISO1012	DR 10ISO1417	DR 10ISO1012	DR 10ISO1417	DRVI 101012N	DR 10ISO1417

## Toshiba



**1**

**FDX2530RPW**  
Wireless



**2**

**FDX3543RPW**  
Wireless

## Flaatz



**3**

**EVS2430W**  
Wireless



**4**

**EVS3643W**  
Wireless

## Carestream



**5**

**DRX-1 Wireless**  
Examion Flexible FPS  
Wireless

## iRay



**6**

**Mars 1417V Wireless**

Producttype	1	2	3	4	5	6
DR-Protectionbox	DR 03ISO1012	DR 03ISO1417	DR 03ISO2430	DR 03ISO1417	DR 03ISO1417	DR 03ISO1417
DR-Skylinebox	DR 05ISO1012	DR 05ISO1417	DR 05ISO2430	DR 05ISO1417	DR 05ISO1417	DR 05ISO1417
DR-Combibox LM	DR 06ISO1012 1	DR 06ISO1417 1	DR 06ISO2430 1	DR 06ISO1417 1	DR 06ISO1417 1	DR 06ISO1417 1
DR-Combibox LM / AP	DR 06ISO1012 2	DR 06ISO1417 2	DR 06ISO2430 2	DR 06ISO1417 2	DR 06ISO1417 2	DR 06ISO1417 2
DR-Combibox SET	DR 06ISO1012 S	DR 06ISO1417 S	DR 06ISO2430 S	DR 06ISO1417 S	DR 06ISO1417 S	DR 06ISO1417 S
DR-Podo-Protectionbox	DR 07ISO1012	DR 07ISO1417	DR 07ISO2430	DR 07ISO1417	DR 07ISO1417	DR 07ISO1417
DR-Cowboybox	DR 08ISO1012	DR 08ISO1417	DR 08ISO2430	DR 08ISO1417	DR 08ISO1417	DR 08ISO1417
DR-Cowboy Package	DR 08ISO1012 P	DR 08ISO1417 P	DR 08ISO2430 P	DR 08ISO1417 P	DR 08ISO1417 P	DR 08ISO1417 P
DR-Sleeve	DR 10ISO1012	DR 10ISO1417	DR 10ISO2430	DR 10ISO1417	DR 10ISO1417	DR 10ISO1417

## Thales



1

Pixium Portable 2430EZ  
Leonardo DR Mini

Wireless

## Idexx



2

Examion DR810  
Idexx Equiview  
Idexx DR 810  
FDXD 810

## Sound Eklin



3

Eklin Mark II Idaho  
Eklin Mark III Idaho



4

Eklin Mark II Texas  
Eklin Mark III Texas

## Varian



5

Paxscan 2520E+  
Leonardo DR2520+  
Protec Rapixx 2520V



6

Paxscan 2520E+  
Cuattro Uno EQ  
Sound Eklin truDR  
Protec Rapixx 2520V

Producttype	1	2	3	4	5	6
DR-Protectionbox	DR 03ISO2430	DRID 03334225	DREK 03IDAHO	DREK 03TEXAS	DRLE 03234528	DRCU 03234528
DR-Skylinebox	DR 05ISO2430	DRID 05334225	DREK 05IDAHO	DREK 05TEXAS	DRLE 05234528	DRCU 05234528
DR-Combibox LM	DR 06ISO2430 1	DRID 06334225 1	DREK 06IDAHO 1	DREK 06TEXAS 1	DRLE 06234528 1	DRCU 06234528 1
DR-Combibox LM / AP	DR 06ISO2430 2	DRID 06334225 2	DREK 06IDAHO 2	DREK 06TEXAS 2	DRLE 06234528 2	DRCU 06234528 2
DR-Combibox SET	DR 06ISO2430 S	DRID 06334225 S	DREK 06IDAHO S	DREK 06TEXAS S	DRLE 06234528 S	DRCU 06234528 S
DR-Podo-Protectionbox	DR 07ISO2430	DRID 07334225	DREK 07IDAHO	DREK 07TEXAS	DRLE 07234528	DRCU 07234528
DR-Cowboybox	DR 08ISO2430	DRID 08334225	DREK 08IDAHO	DREK 08TEXAS	DRLE 08234528	DRCU 08234528
DR-Cowboy Package	DR 08ISO2430 P	DRID 08334225 P	DREK 08IDAHO P	DREK 08TEXAS P	DRLE 08234528 P	DRCU 08234528 P
DR-Sleeve	DR 10ISO2430	DRID 10334225	DREK 10IDAHO	DREK 10TEXAS	DRLE 10234528	DRCU 10234528



# 6

## Ultrasound Solutions



It's not just X-ray accessories Podoblock has to offer. In the realm of ultrasound we also have become a key player. As ultrasound machines have developed at a rapid pace, and thus have become compact and mobile, we offer a gamut of support products.

Take for example the SCANdroid, that has proven itself to be an indispensable and must-have product, also outside the ultrasound sphere.

On the following pages you will find all scanning solutions you'll ever need.

## Cooperflex

A light-weight, sturdy and highly mobile hoof support with a pivot that allows you to change the angle of the foot pedal

Comes with an extension-pipe that gives you more working-height

- Material Base plate: High-density Polyethylene
- Material Foot stand: Aluminium
- Weight: 4,5 kg
- Diameter Base: Ø 40 cm
- Low position: 20 cm
- High position: 40 cm



The extension pipe also serves as a carrying device



## SCAN Wedge

Straight-forward wedge for extra extension to be placed under the hoof

- Material : Rubber 95° Shore
- Dimensions: 180 x 150 x 65-0 mm
- Angle: 20°

## U-scan

Our own standoff, adaptable for a wide range of linear-array transducers

- Thickness sides: 7 mm
- Thickness bottom: 12 mm
- Comes in a protective box



Product	Ordernumber	Price
Cooperflex	CF 10FLSTND	€ 750,-
SCAN Wedge	SW 1050065-0	€ 90,-
U-scan	LAS 10804010	€ 235,-

# SCANdroid

The ideal, height-adjustable scan stool for the mobile veterinarian to transport to the job site or for clinic use

The SCANdroid has great stability with a total of 5 castor wheels, 2 of which with a locking device

The gas spring allows for stageless height-adjustment

- Diameter: 40 cm
- Transport height: 23 cm
- Weight: 9,3 kg



## Transport position



## Position 1: 23 cm

Scanning the foot and lower parts of the leg



## Stageless height adjustment: 38-51 cm

Scanning the front knee and middle parts of the leg



Scanning the stifle and higher parts of the leg



Stow-away room for gas spring and any of your extension(s)

If you require extra height you can order 2 different sizes of optional gas spring extensions with stageless height-adjustment:

- Short: 52 up to 64 cm (stageless)
- Long: 63 up to 75 cm (stageless)



Product	Ordernumber	Price
SCANdroid	SK 1039270	€ 450,-
Extension Short	SK 10EXT150	€ 90,-
Extension Long	SK 10EXT307	€ 90,-

## The SCANster

The optimal ultrasound solution for creating the ideal working height at all times

Especially designed for the portable ultrasound

- Material: Light-weight durable aluminium
- Transport size: 58 x 41 x 26 cm
- Weight: 13,1 kg
- Plateau size: 55 x 37,5 cm
- Maximum load: 15 kg
- 4 locking Ø 80 mm castor wheels



**Position 1: 25 cm**  
For imaging the foot and pastern



**Position 2: 42 cm**  
For imaging the forelimbs and hind limb



**Position 3: 77 cm**  
For imaging the stifle, shoulder and elbow



**Position 4: 112 cm**  
For imaging the hip, sacrum, back, neck and endo-rectal images



With optional proberack



Transport position



Optional with professional probe rack to store a total of 2 ultrasound probes, gel bottle and Standoff PAD

Product	Ordernumber	Price
SCANster	VS 1045321910	€ 975,-
Proberack	CUPH 102P1B1S	€ 100,-



## NewCart

Durable, in height adjustable cart for all small and big machines or equipment, accommodated with gas spring for variable, stageless height-adjustment and storage space in bottom reservoirs

Easy transport position with the pull-out handle, easy storage in vertical position

- Transport size: 42 x 62 x 30 cm
- Weight: 17,2 kg
- Plateau size (inside rim): 57 x 37 cm
- Height-adjustment from 30 to 106 cm
- Maximum load: 15 kg
- 4 swivel, lockable Ø 100 mm castor wheels



**Transport position**



**Position 1:**  
Height 30 cm



**Stageless height adjustment:**  
from 30-106 cm



**Optional ProbeRack**



The Proberack can be attached to any side of the table, at any spot and fits snugly in the bottom casing. It can hold 2 ultrasound probes and a gel bottle

Product	Ordernumber	Price
NewCart	VS 10NCART	€ 1.490,-
Proberack	CUPH 10NCART	€ 100,-

# 7

## Generator Solutions



With our product line of generator gear, you will have no trouble finding a solution that will fit your daily needs and demands. Our Stat-X, JiXtr or X-ray switches will be a true asset to your equine practice, thanks to the all-purpose character of our tools. We can offer standard solutions for every mobile generator.

With regard to your safety, ease-of-use, quality and convenience, you, your staff and your horses will be in good hands with Podoblock equipment!

## TruXRAY™ Switch

Easy prep-and-expose switch for any portable X-ray source available on the market

- Plug-and-play device
- Comes with adapter for your source



Switch size	Ordernumber	Price
152 mm	TXS 10G-2432	€ 250,-
165 mm	TXS 10G-2433	€ 250,-
178 mm	TXS 10G-2434	€ 250,-
190 mm	TXS 10G-2435	€ 250,-
203 mm	TXS 10G-2436	€ 250,-
216 mm	TXS 10G-2437	€ 250,-

**Note:** Please measure the width of your generator grip and quote your generator brand when ordering



Footswitch



Handswitch

## TruXRAY™ Wireless Switch

Wireless foot- or handswitch designed to give the user increased opportunity to get the shot right.

- Does not conflict with any other wireless product
- Requires no pairing
- Plug-and-play device
- Only available in black

Product	Ordernumber	Price
Wireless Hand - Poskom	TXS 10G-2420	€ 495,-
Wireless Hand - EcoRay	TXS 10G-2421	€ 495,-
Wireless Hand - MinXray/Giether	TXS 10G-2422	€ 495,-
Wireless Hand - Porta	TXS 10G-2430	€ 495,-
Wireless Foot - Poskom	TXS 10G-2425	€ 575,-
Wireless Foot - EcoRay	TXS 10G-2426	€ 575,-
Wireless Foot - MinXray/Giether	TXS 10G-2427	€ 575,-
Wireless Foot - Porta	TXS 10G-2428	€ 575,-

**Note:** In case you have a different brand of X-ray source, please contact us. In virtually all cases the switch will fit your source.

## JiXtr

With the JiXtr your X-ray source becomes an easy-to-handle, versatile and highly portable instrument that will cover every angle, spot or height you desire, while working quickly and safely and sparing your back in the process

The JiXtr can be adjusted without pre-defined increments and has a variable height-adjustment. It is compatible with all X-ray devices on the market, with a maximum of 15 kg

- Material: Light-weight aluminium
- Transport size: 12 x 15 x 75 cm
- Low position: 27 x 29 x 69 cm
- High position: 27 x 29 x 201 cm
- Stroke range: 165 cm
- Weight: 5,0 kg
- 'Fix-and-release' system



The rotation of the source in combination with the unlimited tilt of the JiXtr secures an accurate position for Skyline views



Product	Ordernumber	Price
JiXtr	JSM 09JIXTR-14	€ 1.500,-



## The Stat-X Family

With our Stat-X you have the most advanced, all-purpose mobile stand that will fully satisfy your needs in your daily working environment. On the following 3 pages we will show you all the features of the several Stat-X options we have to offer, so you can explore, combine and order the stand with the options of your choice

To make the ordering process easier for you, we've outlined the choice of options and all the basic available models on the following pages, so you can easier decide whether you need a small or a large version of the Stat-X mobile stand. You can then continue to order the extras you desire in order to optimally service your needs for your working environment

If you need further assistance, don't hesitate to call us, so we can inform you more about all the options the Stat-X family has to offer



**Vaquero**



**Espléndido**



**Bucky**



**Infinity**

- 3 swing arm adapter positions
- Swivel range arm: 5° left / right
- All cassette holders are 3D adjustable
- 4 castor wheels



The Stat-X Espléndido is fitted with a new, revolutionary 3D-head

Your generator can be secured in fixed 12° positions with a maximum of 35° left and right



All Stat-X stands are easily folded up in a matter of seconds

It can be transported with the front grip handle on the main tube, while the castor wheels will facilitate moving around

# The Stat-X options

## Swingarm



Available in 2 sizes:

- Swingarm 750 mm Vertical travel 670 mm  
Load capacity: 6,0 - 20,0 kg
- Swingarm 580 mm Vertical travel 500 mm  
Load capacity: 20,1 - 25,0 kg

The quick-release system will rapidly position the swingarm in all desired heights

## Cassette Holder



The cassette holder will be provided, depending on whether you work with an analogue CR- or a DR-system

The cassette holder and the 3D swivel-joint are included

Optional is a tailor-made Grid Compartment

## 3D-Head



The 3D head gives you twelve fixed 12° positions with a maximum of 35° to either side and can be turned totally upside down

The 360° rotation with staged locking of the source secures an accurate position for both skyline views and head X-rays

The 3D head will provide for all portable generators available on the market

## Short legs



All Stat-X stands (except for the Stat-X Espléndido) are available for stationary use with short legs and extra counterweight to create more working space between horse and stand

A counterweight at the back of the stand will prevent it from tipping over

## Color



All Stat-X stands come in one standard color:

- PMS260C Satin



All Stat-X models are CE-marked Class 1 medical products  
European Medical Devices

MD 93/42/EEG  
NL-CA002-2012-24232

# Ordering the Stat-X

## Stat-X Vaquero

- Weight: 22 kg
- Swing arm length: 750 mm
- Vertical travel arm: 700 mm
- Vertical travel stand: 700 mm
- Max height bracket: 209 cm
- Weight portable X-ray source: 6.0 - 20.0 kg



- Swing arm length: 580 mm
- Vertical travel arm: 500 mm
- Vertical travel stand: 700 mm
- Max height bracket: 198 cm
- Weight portable X-ray source: 20.1 - 25.0 kg

- Transport size: L 147 x W 30 x H 32 cm

Optional with short legs (50 cm instead of 90 cm) and extra counterweight (6.5 kg)

Weight X-ray source	Ordernumber	Price
6.0 - 20.0 kg	RSM 09750LW	€ 3.050,-
20.1 - 25.0 kg	RSM 09580HW	€ 3.050,-
Short legs /extra counterweight	RSM 09SLEW	€ 235,-

## Stat-X Bucky

Includes a cassette holder for analogue, CR or DR systems

- Weight: 22 kg
- Weight X-ray cassette / Panel: 4.0 - 20.0 kg
- Swing arm length: 750 mm
- Vertical travel arm: 700 mm
- Vertical travel stand: 410 mm
- Max height bracket: 185 cm



- Transport size: L 110 x W 30 x H 32 cm

Optional with short legs (50 cm instead of 90 cm) and extra counterweight (6.5 kg)

Weight X-ray cass/panel	Ordernumber	Price
4.0 - 20.0 kg	BSM 09750LW	€ 3.550,-
Short legs /extra counterweight	RSM 09SLEW	€ 235,-

## Stat-X Espléndido

- Weight: 22 kg
- Weight portable X-ray source: 6.0 - 20.0 kg
- Swing arm length: 750 mm
- Vertical travel arm: 700 mm
- Vertical travel stand: 410 mm
- Max. height bracket: 195 cm



- Transport size: L 114 x W 30 x H 32 cm

Weight X-ray source	Ordernumber	Price
6.0 - 20.0 kg	ESM 09750LW	€ 3.750,-

## Stat-X Infinity

Includes a cassette holder for analogue, CR or DR systems

Technical features comparable to Stat-X Vaquero

- Weight: 26 kg
- Weight portable X-ray source: 6.0 - 20.0 kg  
20.1 - 25.0 kg
- Transport size: L 147 x W 30 x H 32 cm



Optional with short legs (50 cm instead of 90 cm) and extra counterweight (6.5 kg)

Weight X-ray source	Ordernumber	Price
6.0 - 20.0 kg	ISM 09750LW	€ 4.400,-
20.1 - 25.0 kg	ISM 09580HW	€ 4.400,-
Short legs /extra counterweight	RSM 09SLEW	€ 235,-

# 8

## Protection & Safety



Radiation protection is pivotal to all X-ray exams. It's the alpha and the omega. Fortunately we offer a complete range of products in order to optimally protect our users against radiation. With personalized sizes and colours we have tried to make our product line suitable for all our customers.

We have solutions for all environments, be they stationary at your clinic or in a mobile environment - on any location really.

## Thyroid Shield

- Lightweight nylon
- Velcro fastener
- Pb 0.5 mm
- One size fits all
- Standard delivered in color #13 from stock
- other colors with extra charge



Product	Ordernumber	Price
Thyroid Shield	RS 10BL05M	€ 80,-

For available colors and prints, see table on page 43

## X-ray Glove

- Flexible vinyl glove
- Pb 0.5 mm



Size	Ordernumber	Price
Small	RH 10LF05S	€ 210,-
Medium	RH 10LF05M	€ 210,-

## X-ray Protection Goggles

- Pb 0.5 mm
- Comes in a protective case



Product	Ordernumber	Price
X-ray Protection Goggles	RB 10OP05	€ 185,-



## X-ray Apron

- Pb 0.50 mm Front
- Pb 0.25 mm Back
- Snap fasteners



Size		Ordernumber	Price
Small	50 x 110 cm	RJ 10CLI05S	€ 190,-
Medium	60 x 120 cm	RJ 10CLI05M	€ 210,-
Large	75 x 120 cm	RJ 10CLI05L	€ 240,-

Standard delivered in color #13 (from stock), lead-free version available at request  
For other available colors and prints see table on opposite page

## X-ray Apron "de Luxe"

With integrated Thyroid shield protection

- Pb 0.50 mm Front
- Pb 0.25 mm Back
- Pb 0.50 mm Thyroid
- Snap fasteners



Size		Ordernumber	Price
Small	50 x 110 cm	RJ 10CLI05STS	€ 260,-
Medium	60 x 120 cm	RJ 10CLI05MTS	€ 275,-
Large	75 x 120 cm	RJ 10CLI05LTS	€ 310,-

Standard delivered in color #13 (from stock), lead-free version available at request  
For other available colors and prints see table on opposite page

## X-ray Apron Vest & Skirt

Reduces stress on your back to a minimum

Equipped with shoulder pads

- Pb 0.50 mm Front
- Pb 0.25 mm Back
- Snap fasteners



### Available sizes:

(Vestlength - Skirtlength - Circumference - Total length - all in cm)

Size		Ordernumber	Price
Small	60-60-107-105	RJ 10VS05S	€ 625,-
Medium	63-65-117-110	RJ 10VS05M	€ 710,-
Large	65-75-127-120	RJ 10VS05L	€ 760,-

Standard delivered in color #13 (from stock), lead-free version available at request  
For other available colors and prints see table on opposite page

## Color card X-ray Apron, X-ray Apron "de Luxe" and X-ray Apron Vest & Skirt (also applies to 'VectorX Light Lead' orders)

Please place the requested color number on your order behind the Apron article number



Nr. 1 - Dark Green



Nr. 13 - Dark Blue



Nr. 4-1 - Dark Purple



Nr. 12 - Leopard



Nr. 6-1 - Dark Red



Nr. 15 Camouflage



Nr. 8 - Dark Grey



Nr. 5-5 - Python



Nr. 9 - Yellow



Nr. 2-1 - Red with Horse



## Warning sign

Mobile warning sign that radiation is being released

Available languages

- English
- German
- Dutch
- Danish
- French



Other languages on request

Language	Ordernumber	Price
English	WS 10600EN	€ 50,-
German	WS 10600DE	€ 50,-
Dutch	WS 10600NL	€ 50,-
Danish	WS 10600DK	€ 50,-
French	WS 10600FR	€ 50,-

## Barrier tape

- Two-sided print
- Easy tear-off
- Width: 75 mm
- 250 m per roll



Product	Ordernumber	Price
1 Roll	WT 1025075-1	€ 50,-
10 Rolls	WT 1025075-10	€ 225,-

## Mobile Lead Shield

A sturdy protective shield made of anodized aluminium, lined with an impact-resistant plastic layer on both sides

Comes with four Ø80 mm castor wheels, two of which with brakes

24 x 30 cm lead glass window

- Height: 190 cm
- Width: 80 cm
- Pb 2.0 mm



Product	Ordernumber	Price
Mobile Lead Shield	MLS 1025076986	€ 2.350,-

Glove Rack

- 3 pair gloves
- Wall mount




Apron Rack Wood

- 1 garment
- Wood



Apron Rack

- 5 garments
- Wall mount

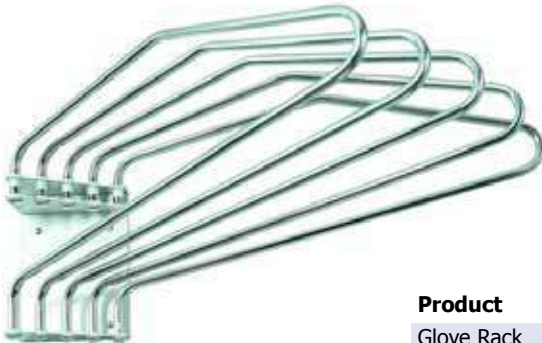

Product	Ordernumber	Price
Glove Rack	RR 10CH03	€ 225,-
Apron Rack Wood	AH 10WMI	€ 12,50
Apron Rack	RR 10CH05	€ 325,-
Apron Rack Round	RR 10CH08	€ 375,-

Apron Rack Round

- 8 garments
- Mobile
- Height: 165 cm



# 9

# X-ray Markers



In spite of the digital revolution the use of markers is still commonplace

Isn't it convenient to be able to carry on doing this in the 21<sup>st</sup> century as if it were the analog era? Sometimes there's just nothing wrong with being 'Old School'



## X-ray Clip

Stainless steel clip for  
analogue and CR-cassettes

- 55 x 21 mm
- Gap: 18 mm



Product	Ordernumber	Price
L / R	RC 10LR18	€ 10,-
LF / RF	RC 10LFRF18	€ 10,-
LV / RV	RC 10LVRV18	€ 10,-
LH / RH	RC 10LHRH18	€ 10,-
LAT / MED	RC 10LATMED18	€ 10,-



## X-ray Clip DR

Multiclip that clenches around  
the edges of our boxes / sleeves

- 95 x 25 mm
- Gap: 27 mm

Product	Ordernumber	Price
L / R	RC 10LR27	€ 10,-
LF / RF	RC 10LFRF27	€ 10,-
LV / RV	RC 10LVRV27	€ 10,-
LH / RH	RC 10LHRH27	€ 10,-
LAT / MED	RC 10LATMED27	€ 10,-

## Multimarker

Stainless steel all-in-one  
marker with suction cup

- Ø 50 mm



Product	Ordernumber	Price
LF / RF / LH / RH	RCM 10ENZN	€ 25,-
LV / RV / LH / RH	RCM 10DEZN	€ 25,-
LAT / MED / DV / VD	RCM 10LMZN	€ 25,-



## Multiclip DR

Multiclip that clenches around  
the edges of our boxes or  
sleeves

- Ø 50 mm
- Gap 27 mm

Product	Ordernumber	Price
LF / RF / LH / RH	RCM 10ENDR	€ 25,-
LV / RV / LH / RH	RCM 10EDDR	€ 25,-
LAT / MED / DV / VD	RCM 10LMDR	€ 25,-

## Sarah marker

Marker with integrated label  
to always be sure about the  
horizontal plane



Product	Ordernumber	Price
L / R with suction cup	SC 10LRZN	€ 60,-
L / R with stainless steel clip	SC 10LRDR	€ 60,-



## Buckymarker

Marker with suction cup or clip  
for your cassette holder to  
provide information about the  
horizontal plane

- 360° rotation
- Size: 48 x 44 x 4 mm

Product	Ordernumber	Price
With suction cup	MKB 1048444	€ 40,-
With stainless steel clip	MKB 1048444DR	€ 40,-

# 10

## Special X-ray Solutions



After our early-day inventing and producing of a number of X-ray accessories, we received a steady flow of demand for specialized and customized products to enhance and professionalize the quality of X-rays.

All these inventions have meanwhile become mainstream in our product line and can usually be delivered from stock.

## Aluminium Grid

Aluminium cover with parallel lines

- 70 lines / cm (178 lines / inch)
- Focus: Infinity
- Ratio: 8:1
- Only to be used with Analogue and CR systems



Product	Ordernumber	Price
18 x 24 Film size: 180 x 240 mm	RL 10182415	POR*
18 x 24 Cass. size: 200 x 260 mm	RL 10202615	POR*
24 x 30 Film size: 240 x 300 mm	RL 10243015	POR*
24 x 30 Cass. size: 260 x 320 mm	RL 10263215	POR*
35 x 43 Film size: 356 x 432 mm	RL 10354315	POR*
35 x 43 Cass. size: 381 x 457 mm	RL 10384615	POR*

POR = Price on request

## Aluminium Wedge

Filter the radiation with high purity Homogenous quality aluminium for better back / neck X-rays

2 small wedges with Velcro for use on your X-ray generator

Large concave wedge by Prof. J.M. Denoix, provided with counterweight for easy fixation in all positions



Product	Ordernumber	Price
195 x 75 mm / 23°	AK 101957530	€ 375,-
195 x 50 mm / 45°	AK 101955045	€ 375,-
400 x 200 x 40-0 mm	AK 1040025040	€ 1.250,-

## Backscatter Lead

Separate foil with sealed lead to avoid backscatter

Fast and convenient application behind any X-ray panel

Fully covered lead Pb 0.15 mm



Size	Ordernumber	Price
268 x 329 mm	BL 10ISO2430	€ 150,-
282 x 333 mm	BL 10ISO1012	€ 150,-
385 x 460 mm	BL 10ISO1417	€ 150,-

Sizes other than above can be ordered at request

## EZ Paste



A Barium Sulphate, radiopaque paste that provides excellent radiographic visualization of the true hoof wall

Product	Ordernumber	Price
EZ Paste (454 gr)	BP 10EZPASTE	€ 115,-

## Aiangle

A great device to measure the exact angle. It is easy to place on both your generator and Cassette Holder.

Analog

- Measuring range: 360°
- Velcro base

Digital

- Measuring range: 360°
- Magnetic base, tilt display
- 2 AAA batteries included



Product	Ordernumber	Price
Analog	AI 10AN1C	€ 45,-
Digital	AI 10RAS2C	€ 98,-



Changing a horseshoe is ideally done by a farrier, but when an equine vet chooses to do so, we have a number of professional and ingenious tools for you, lined out in this chapter.

You will also find for instance our Hoofplaster, that is already used in many practices, in order to ameliorate navicular bone X-rays.

## VET Tools

A comprehensive set of tools to get you going, quickly and safely



### VET Tools Set includes:

- Apron "VET"
- Soft-faced Hammer
- Clench Cutter
- Straight Blade
- Loop Knife
- Nail Puller
- Pull-off
- Adjustable Hooftester
- Canvas Storage Bag

#### Product

Apron "VET"  
VET Tools Set

#### Ordernumber

VV 10VET  
VV 10SET



## Adjustable Hooftester

Easy-to-use, adjustable hooftester

- Aluminum

#### Product

Adjustable Hooftester

#### Ordernumber

VV 10WERK07



## 123 Shoeing Kit

Use this indispensable high-quality farrier-hammer with the clinch cutter for removing 'hidden' nails during de-shoeing

- Shoe Puller
- Hammer
- Clinch Cutter

#### Product

123 Shoeing Kit

#### Ordernumber

PZT 10123SK

## Sørender

For use directly between scissor teeth in combination with full anaesthetic with a spare Ø 1.18" or Ø 1.97" hole for tube transit

- Material: PP
- Biothane strap



#### Size

M: 3.94" x 3.94" x 1.97" / Ø1.18  
L: 4.72" x 4.53" x 2.17" / Ø1.97

#### Ordernumber

AB 1010030  
AB 1012050



## Hoofplaster

For shadow-free hoof recordings

#### Content

5 x 8.8 oz  
25 x 8.8 oz  
50 x 8.8 oz

#### Ordernumber

HP 10BL25005  
HP 10BL25025  
HP 10BL25050





## Podoblock B.V.

Verlengde Nijverheidsweg 10  
9482 WT Tynaarlo  
The Netherlands

Office 0031-508 200 257  
Fax 0031-508 200 212  
Alex Berends 0031-651 828 117

i: [www.podoblock.com](http://www.podoblock.com)  
e: [info@podoblock.com](mailto:info@podoblock.com)